

# Reference Guide for Radiologists MRI INFORMATION CHECKLIST



# 1991 Clarion C1.0 Clarion C1.2 1995 2000 Clarion CII HiRes 90K 2003 2013 HiRes 90K Advantage HiRes Ultra 2016 2018 HiRes Ultra 3D

## Reference Guide for Radiologists

#### MRI INFORMATION CHECKLIST FOR ADVANCED BIONICS COCHLEAR IMPLANTS PATIENTS

For Europe Only

The Instructions in the present document only apply to countries with CE Mark approval.

A summary of the MRI Safety Information for all AB Cochlear Implants is provided below. Please proceed to the checklist for the applicable implant for detailed information.

IMPLANT TYPE	MRI FIELD STRENGTH (T)	SPATIAL GRADIENT FIELD (T/m)		MAX HEAD SAR (W/kg)	MAX BODY SAR (W/kg)
C1.0*	MRI contraindicated	no	ne	none	none
C1.2*	MRI contraindicated	no	ne	none	none
CII*	MRI contraindicated	none		none	none
HiRes 90K*	1.5T	2.5 T/m		≤1.0 W/kg	≤1.7 W/kg
HiRes 90K Advantage	1.5T	2.5	2.5 T/m		≤1.7 W/kg
HiRes Ultra	1.5T	3.47 T/m	13.90 T/m***	≤ 3.2 W/kg	≤ 2.0 W/kg
nines Oitra	3.0T**	6.9 T/m		≤ 2.6 W/kg	≤ 2.0 W/kg
LliDag Llitra 2D	1.5T	20 T/M		≤ 3.2 W/kg	≤ 2.0 W/kg
HiRes Ultra 3D	3.0T	20 T/M		≤ 2.6 W/kg	≤ 2.0 W/kg

<sup>\*</sup> This device is no longer sold in the EU or North America

For additional information regarding the use of an MRI scanner with a HiRes Ultra 3D, please contact Advanced Bionics Technical Support at technicalservices@advancedbionics.com or visit www.advancedbionics.com/mri.

<sup>\*\*</sup> For MRI the magnet has to be removed surgically

<sup>\*\*\*</sup> With magnet removed

### HiRes Ultra 3D



The HiRes™ Ultra 3D cochlear implant has an innovative designed multi-magnet assembly that allows safe MRI scanning with the multi-magnet assembly in place, without angular restrictions of the head, or a bandaging protocol.

#### MRI Checklist

#### Step 1: Pre-Scan Patient Screening

Before scheduling an MRI, it's vital to contact the patient's CI clinic or surgeon to answer these questions and fully understand their situation and scanning need<sup>1</sup>.

	A. Is the user bilaterally implanted (has implants on both left and right ears) with HiRes Ultra 3D			
	☐ YES	Proceed to question B.		
	□ NO	The patient has only one cochlear implant and it is HiRes Ultra 3D. Proceed to question B.		
	□ NO	The patient is bilaterally implanted but has HiRes Ultra 3D on one side and a different CI on the other.		
		Verify that both implants are conditionally approved for MRI before proceeding with a scan <sup>2</sup> . <b>NOTE</b> : The rest of this checklist will not be applicable to the case and it's important to follow the necessary programming, bandaging, surgery and angle protocols for the other CI to ensure patient safety.		
		e patient's doctor requesting a reduced image-artifact size around the cochlear implant site ample, some head, brain or neck scans <sup>3</sup> )?		
	□ NO	Proceed to Step 2: Prep Patient for Scan.		
	☐ YES	Prior to MRI Scan:		
		☐ Order a Temporary Non-magnetic Plug (CI-1420) for every HiRes Ultra 3D implant the patient has (2 for bilateral implants, 1 for a single implant).		
		☐ Order a HiRes Ultra 3D Replacement multi-magnet assembly (CI-1419) for every HiRes Ultra 3D implant the patient has (2 for bilateral implants, 1 for a single implant).		
		☐ Schedule surgery with the patient's cochlear implant surgeon to remove the multi-magnet assembly <sup>4</sup> and possibly replace it with the Temporary Non-Magnetic Plug(s) (CI-1420).		
		☐ After the scan, schedule surgery with the patient's cochlear implant surgeon to replace the Temporary Non-Magnetic Plug (CI-1420) <sup>4</sup> with the HiRes Ultra 3D Replacement multi-magnet assembly (CI-1419).		
S	tep 2: Pa	atient Preparation on Day of Scanning		
		are that it has at least been 2-4 weeks since the patient cochlear implant surgery. <b>NOTE:</b> This is the recommended mum duration in order to allow any inflammation from the cochlear implant to subside prior to MRI scan.		
	☐ Verify	y the patient does not have any fever.		
	☐ Cour	nsel the patient regarding the possibility of auditory sensations during the scan.		

Remove patient's sound processor and headpiece before entering a room containing an MRI scanner. NOTE: The

external sound processor and headpiece are MR unsafe.

Proceed with MRI scan in adherence with the following important scanning guidelines:

#### MRI Scan Safety Guidelines

Testing has demonstrated that the HiRes Ultra 3D cochlear implant is MR Conditional. Unilateral and bilateral recipients with this device with multi-magnet assembly in place can be safely scanned in a horizontal closed-bore quadrature coil MR system meeting the following conditions:

MRI FIELD STRENGTH	1.5 T	3.0 T
Maximum Spatial Field Gradient	20 T/m	20 T/m
RMS Gradient Field	34.4 T/s	34.4 T/s
Peak Slew Rate	200 T/m/s	200 T/m/s
Maximum whole body averaged SAR	2.0 W/kg	2.0 W/kg
Maximum head averaged SAR	3.2 W/kg	2.6 W/kg

When tested under the scan conditions defined above the HiRes Ultra 3D cochlear implant produced a maximum temperature rise of <3°C after 15 minutes of continuous scanning.

**NOTE**: During the scan, patients might perceive auditory sensations. Adequate counseling of the patient is advised prior to performing the MRI. The likelihood and intensity of the auditory sensations can be reduced by selecting sequences with a lower Specific Absorption Rate (SAR) and slower gradient slew rates.

The largest artifact at 3.0T is >15cm when imaged under spin echo and gradient echo sequences. Smaller artifacts are possible using different scanning parameters only.

#### Step 3: Post-Scan Patient Management

#### Was the multi-magnet assembly removed prior to MRI scan?

☐ YES Confirm multi-magnet assembly replacement surgery schedule with the cochlear implant surgeon and patient.

The cochlear implant user may reattach and start using the external components (sound processor and headpiece) upon exiting the MRI field zone.

For additional information regarding the use of an MRI scanner with a HiRes Ultra 3D cochlear implant, please contact Advanced Bionics Technical Support at <a href="mailto:technicalservices@advancedbionics.com">technicalservices@advancedbionics.com</a> or visit <a href="mailto:www.advancedbionics.com">www.advancedbionics.com</a>/mri.

5

### HiRes Ultra



The HiRes™ Ultra cochlear implant is indicated for 3.0 Tesla MRI with magnet removal and 1.5 Tesla MRI scanning with the magnet in and head bandaging protocol.

#### MRI Checklist

#### Step 1: Pre-Scan Patient Screening

Before scheduling an MRI, it's vital to contact the patient's CI clinic or surgeon to answer these questions and fully understand their situation and scanning need.

A. Is the user bilaterally implanted with HiRes Ultra?		
☐ YES Proceed to question B <sup>1,2</sup> .		
□ NO	The patient has only one cochlear implant and it is HiRes Ultra. Proceed to question B <sup>1,2</sup> .	
□NO	The patient is bilaterally implanted and has HiRes Ultra only on one side.	
	Verify that both implants are conditionally approved for MRI before proceeding with a scan <sup>1</sup> . <b>NOTE:</b> The rest of this checklist will not be applicable to the case and it's important to follow the necessary programming, bandaging, surgery and angle protocols for the other non-HiRes Ultra CI to ensure patient safety.	
B. Is the p	atient's doctor requesting a 1.5 Tesla MRI?	
☐ YES	Proceed to head bandaging protocol on pages 10 and 11.	
C. Is the patient's doctor requesting a reduced image artifact size around the cochlear implant (for example, head and neck scans <sup>3</sup> )?		
□ NO	Proceed to head bandaging protocol on pages 10 and 11.	
	patient's doctor requesting a 3.0 Tesla MRI or a reduced image artifact size around the mplant site (for example, head and neck scans <sup>3</sup> )?	
☐ YES	STOP, magnet removal is needed. Please observe the following steps.	
	Prior to MRI Scan:	
<ul> <li>Order a Temporary Non-magnetic Plug (CI-1411) for every HiRes Ultra implant the patient h bilateral implants, 1 for a single implant).</li> </ul>		
	Order a HiRes Ultra Replacement magnet (CI-1413) for every HiRes Ultra implant the patient has (2 for bilateral implants, 1 for a single implant).	
	Schedule surgery with the patient's cochlear implant surgeon to remove the magnet(s) and possibly replace it with the Temporary Non-Magnetic Plug(s) (CI-1411).	
	□ YES □ NO □ NO □ NO □ NO □ YES C. Is the p (for examp □ NO D. Is the p cochlear in □ YES	

After the scan, schedule surgery with the patient's cochlear implant surgeon to replace the Temporary

Non-Magnetic Plug(s) with the HiRes Ultra Replacement magnet (CI-1413).

The measured range of MRI device image artifact radius extending from the HiRes Ultra cochlear implant for spin and gradient echo sequences and all imaging planes are provided in tables A and B.

Table A: Artifact at 3.0T MRI Field Strength (magnet removed - Spin echo sequence)

IMPLANTED	MAGNET	ARTIFACT RANGE
Unilaterally	Removed	4.7 cm
Bilaterally	Removed	4.7 cm

Table B: Artifact at 1.5T MRI Field Strength (Gradient echo sequence)

IMPLANTED	MAGNET	ARTIFACT RANGE
Unilotorally	In place	~ 7.9 cm
Unilaterally	Removed	~ 3.1 cm
Dilatorally	In place	> 9.5 cm
Bilaterally	Removed	~ 4.2 cm

#### Step 2: Patient Preparation on Day of Scanning

Ensure that	t it has at least been 2-4	4 weeks since the patie	nt cochlear implant surge	ry. NOTE: This is the reco	mmended
minimum (	duration in order to allo	ow any inflammation fro	om the cochlear implant t	o subside prior to MRI s	can.

- Verify the patient does not have any fever.
- Counsel the patient regarding the possibility of auditory sensations during the scan.
- Remove patient's sound processor and headpiece before entering a room containing an MRI scanner. **NOTE**: The external sound processor and headpiece are MR unsafe.

#### Important MRI Scan Safety Guidelines

Testing has demonstrated that the HiRes Ultra cochlear implant is MR Conditional. Unilateral and bilateral recipients with this device can be safely scanned in a horizontal closed-bore quadrature coil MR system meeting the following conditions:

MRI FIELD STRENGTH	1.5 T	3.0 T
Maximum Spatial Field Gradient	3.47 T/m (13.90 T/m with magnet removed)	6.9 T/m (with magnet removed)
RMS Gradient Field	34.4 T/s	34.4 T/s
Peak Slew Rate	200 T/m/s	200 T/m/s
Maximum whole body averaged SAR	2.0 W/kg	2.0 W/kg
Maximum head averaged SAR	3.2 W/kg	2.6 W/kg

When tested under the scan conditions defined above the HiRes Ultra cochlear implant produced a maximum temperature rise of <3°C after 15 minutes of continuous scanning and continuous 3.0T scanning.

**NOTE**: During the scan, pateints might perceive auditory sensations. Adequate counseling of the patient is advised prior to performing the MRI. The likelihood and intensity of the auditory sensations can be reduced by selecting sequences with a lower Specific Absorption Rate (SAR) and slower gradient slew rates.

#### Step 3: Post-Scan Patient Management

#### Was the magnet removed prior to MRI scan?

- YES Confirm magnet replacement surgery schedule with the cochlear implant surgeon and patient.
- NO The cochlear implant user may reattach and start using the external components (sound processor and headpiece) upon exiting the MRI field zone.

6

### HiRes 90K Family

for bilateral implants, 1 for a single implant).

replace it with the Temporary Non-Magnetic Plug(s) (CI-1410).

Non-Magnetic Plug(s) with the HiRes 90K Replacement magnet (CI-1412).



The HiRes 90K<sup>™</sup> family devices are indicated for 1.5 Tesla MRI scanning with the magnet in place and head bandaging protocol, or with magnet removed.

#### MRI Checklist

#### Step 1: Pre-Scan Patient Screening

Before scheduling an MRI, it's vital to contact the patient's CI clinic or surgeon to answer these questions and fully understand their situation and scanning need.

nderstand the	eir situation and scanning need.
A. Is the u	user bilaterally implanted with cochlear implants of the HiRes 90K family?
☐ YES	Proceed to question B <sup>1,2</sup> .
□NO	The patient has only one cochlear implant and it is one of the HiRes 90K family. Proceed to question B1.2.
□ NO	The patient is bilaterally implanted and has one of the HiRes 90K family devices on one side.
	Verify that both implants are conditionally approved for MRI before proceeding with a scan <sup>1</sup> . <b>NOTE</b> : The rest of this checklist will not be applicable to the case and it's important to follow the necessary programming, bandaging, surgery and angle protocols for the other non-HiRes 90K CI to ensure patient safety.
B. Is the p	patient's doctor requesting a 1.5 Tesla MRI?
☐ YES	Proceed to head bandaging protocol on pages 10 and 11.
	patient's doctor requesting a reduced image artifact size around the cochlear implant site ple, head and neck scans <sup>3</sup> )?
□ NO	Proceed to head bandaging protocol on pages 10 and 11.
☐ YES	STOP, magnet removal is needed. Please observe the following steps.
	Prior to MRI Scan:
	Order a Temporary Non-magnetic Plug (CI-1410) for every HiRes 90K family device the patient has (2 for bilateral implants, 1 for a single implant).
	Order a HiRes 90K Replacement magnet (CI-1412) for every HiRes 90K family device the patient has (2

Schedule surgery with the patient's cochlear implant surgeon to remove the magnet(s) and possibly

After the scan, schedule surgery with the patient's cochlear implant surgeon to replace the Temporary

The measured range of MRI device image artifact radius extending from the HiRes 90K Advantage cochlear implant for a gradient echo sequence and all imaging planes are provided in table C.

Table C: Artifact at 1.5T MRI Field Strength

IMPLANTED	MAGNET	ARTIFACT RANGE
Unilaterally	In place	~ 11.5 cm
	Removed*	~ 6 cm

<sup>\*</sup> With temporary non-magnetic plug

These artifacts may result in a loss of diagnostic information in the implant vicinity.

#### Step 2: Patient Preparation on Day of Scanning

Ensure that it has at least been 2-4 weeks since the patient cochlear implant surgery. <b>NOTE</b> : This is the recommended
minimum duration in order to allow any inflammation from the cochlear implant to subside prior to MRI scan.

☐ Verify the patient does not have any fever.

Counsel the patient regarding the possibility of auditory sensations, pain, pressure or discomfort during the scan.

Remove patient's sound processor and headpiece before entering a room containing an MRI scanner. **NOTE**: The external sound processor and headpiece are MR unsafe.

#### Important MRI Scan Safety Guidelines

Testing has demonstrated that the devices of the HiRes 90K family are MR Conditional. Unilateral and bilateral recipients with these devices can be safely scanned in a horizontal closed-bore quadrature coil MR system meeting the following conditions:

MRI FIELD STRENGTH	1.5 T	3.0 T
Maximum Spatial Field Gradient	2.5 T/m*	none
RMS Gradient Field	30 T/s	none
Peak Slew Rate	150 T/m/s	none
Maximum whole body averaged SAR	≤ 1.7 W/kg	none
Maximum head averaged SAR	≤ 1.0 W/kg	none

<sup>\*</sup> With temporary non-magnetic plug

When tested under scan conditions defined above with the patient's head in the center of the MR scanner, the HiRes 90K Advantage implant produced a maximum temperature rise of <2.4°C after 15 minutes of continuous scanning when using a body coil and whole body averaged SAR of 1.7 W/kg. NOTE: During the scan, patients might perceive auditory sensations. Adequate counseling of the patient is advised prior to performing the MRI. The likelihood and intensity of the auditory sensations can be reduced by selecting sequences with a lower Specific Absorption Rate (SAR) and slower gradient slew rates.

#### Step 3: Post-Scan Patient Management

#### Was the magnet removed prior to MRI scan?

☐ YES Confirm magnet replacement surgery schedule with the cochlear implant surgeon and patient.

NO The cochlear implant user may reattach and start using the external components (sound processor and headpiece) upon exiting the MRI field zone.

## Head Bandaging Protocol

#### **Required Material**

- Hurt-Free Wrap Hospital Grade, Latex Free, 5 cm width, Johnson & Johnson (or equivalent)
- Flexible tape, Latex Free, 2 2.5 cm width, 3M Nexcare (or equivalent)
- Patient Headpiece
- MRI Antenna Coil Cover (CI-7521) Advanced Bionics
- Marking Pen

The MRI Antenna Coil Cover is intended to be used along with a bandaging protocol to allow the patient to undergo an MRI procedure with the magnet left in place. Instructions for use must be followed, including the recommended bandaging protocol.

#### Compression Dressing prior to MRI Procedure with HiRes 90K Users

This document will provide general guidance for applying appropriate compression dressing over the site of the implanted device (and magnet) in order to perform an MRI procedure.

1. Place patient in sitting position to allow access to the implant site. **NOTE:** If the patient is bilaterally implanted, prepare to locate both implant devices using the following steps. Both implant locations must have compression dressing placed over their magnets.

#### Secure MRI Antenna Coil Cover over implant magnet site

- 2. Place the patient headpiece, with cable removed, over the implant site (Fig. A). The magnets will hold the headpiece in place.
- 3. Cut a piece of the Coach bandage that is long enough to wrap around the head once.
- 4. Wrap tight this piece around the head, so that the bandage covers the patient headpiece. (Fig. B)
- 5. Outline the position of the headpiece on the bandage using a marker or pen. (Fig. C)
- 6. Slip out the patient headpiece, but keep the bandage in place. (Fig. D)
- 7. Slide the MRI Antenna Coil Cover under the bandage (Fig. E), lining it up with the outline of the headpiece. (Fig. F)



A. Patient heapiece magnetically attached



D. Removing headpiece with bandage in place



B. Patient heapiece under bandage



E. Inserting MRI Antenna Coil Cover



C. Marked bandage showing outline of headpiece




F. MRI Antenna Coil Cover under marked outline

#### Measure head size and bandage length needed for compression wrapping


- 8. Take the remaining Coach bandage roll and wrap the bandage around the head once, without stretching.
- 9. Mark the location on the bandage that is one full wrap around the head. This is the head circumference. (Fig. G)
- 10. Unwrap the piece of bandage with the marked head circumference, and place it on a flat surface.
- 11. Unroll the remaining bandage roll.
- 12. Fold over the bandage start and crease at the marked head circumference. (Fig. H)
- 13. Cut or tear the remaining bandage where it overlaps with the bandage start. (Fig. I) The resulting bandage piece length is twice the head circumference.



G. Remaining Coach bandage roll wrapped once around head

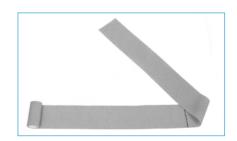

H. Marked head circumference, Fold and crease bandage

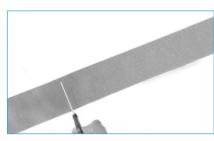

I. Remaining bandage cut where it overlaps with the bandage start

#### Apply compression bandaging

- 14. Wrap this cut piece, this time very tightly, by stretching the marked line an additional half turn around the head. (Fig. J) This ensures 150% extension of the bandage. **NOTE**: After 1.5 turns, the marked line should be on the opposite side of the head from where you started the wrap.
- 15. Continue wrapping the bandage at 150% extension for an additional 1.5 turns, resulting in 3 full turns total. (Fig. K) **NOTE:** After 3 turns, the bandage should end where you started the wrap (further is acceptable also).



J. Tight wrapping, stretching the marked line an additional half turn



K. At least 3 full turns to be competed

#### After the MRI

Once the MRI procedure has been completed, gently remove the coach wrap, tape, and the MRI antenna coil cover. Discard all components including the MRI antenna coil cover.

Please go back to step 2.

10 11



#### Advanced Bionics AG

Laubisrütistrasse 28, 8712 Stäfa, Switzerland T: +41.58.928.78.00 F: +41.58.928.78.90 info.switzerland@AdvancedBionics.com

#### Advanced Bionics LLC

28515 Westinghouse Place Valencia, CA 91355, United States T: +1.877.829.0026 T: +1.661.362.1400 F: +1.661.362.1500 info.us@AdvancedBionics.com

For information on additional AB locations, please visit AdvancedBionics.com/contact

AB - A Sonova brand

<sup>1.</sup> Failure to do so can lead to device movement, device damage, multi-magnet assembly movement, patient discomfort, or trauma and pain to the patient.

If cochlear implant models on either side are different, the MRI safety criteria that is most restrictive of the two implant models must be applied at the discretion of a qualified MRI professional. MRI procedures are contraindicated for CLARION (C1 and CII) cochlear implant recipients. For more information please contact Advanced Bionics Technical Support at <u>technicalservices@advancedbionics.com</u> or visit <u>www.advancedbionics.com/mri</u>.

<sup>3.</sup> Reduced-image artifact scans may require that the implanted multi-magnet assembly be removed prior to scanning.

<sup>4.</sup> Instructions for removal and replacement are provided with CI-1419 and CI-1420.